

vyvo™

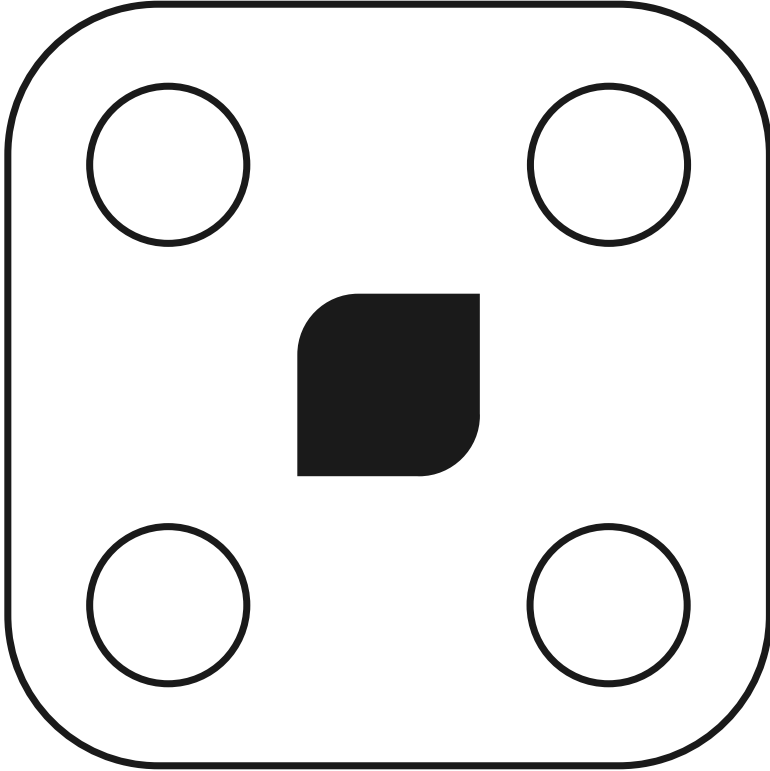leggera

Tracks Weight, BMI, BMR and more using BioElectrical Impedance.
All in just one step.

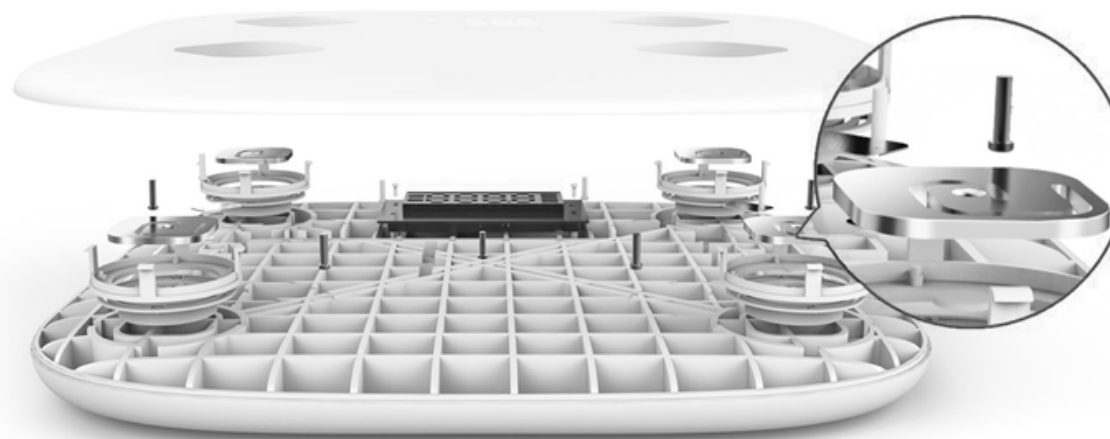QUICK START GUIDE

VYVO LEGGERA scale utilizes BIA (Bio-electrical Impedance Analysis) technology.

Bio-electrical impedance technology makes the measurements more accurate by taking advantage of the conductive composition of muscle and water contents within the body.

The device sends harmless signals through the body to measure the overall body-impedance which is increased due to fat composition and lowered by the conductive muscles and water in the body.


1

1. PACKING LIST


The box contains the following items:

Smart Bluetooth Scale

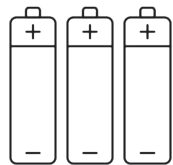Batteries: 3*1.5V AAA

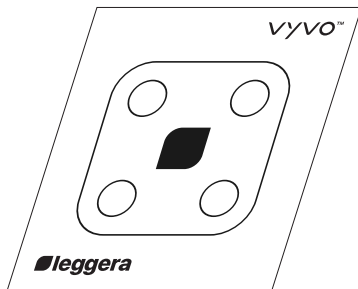Quick Start Manual

2

2. HOW TO OPERATE

2.1. APP Download

Download the VYVO Smart App from the App Store, Google Play Store or VYVOAppSTORE.

Note: To use "VYVO Smart App" on an iPad, you need to select "iPhone only" when downloading it from the App Store.

2.2. Compatibility & Support

Support Bluetooth 4.0.
Android 4.3 or above, iOS 8.0 or above.

3

3. SPECIFICATION PARAMETER

2.3. Connection

1. After downloading the App, please turn on Location Service, enable GPS permission of your mobile phone prior to turn on Bluetooth connection.

2. Register your own ID, then input your information (user name, gender, preferred weight unit, birthday, height).

2.4. View Data

1. After the scale display show three times and stop to your body weight result. That is mean the body composition analysis are measured, please click "Yes" to receive them to your app. You can check the standard of your each body composition by tapping on the data.

2. Click into the history, you can see your history data in graph orst.

| | |
|----------------------|--|
| Product model | LEGGERA |
| Size | 260 × 260 × 27 mm |
| Capacity | 6.6-330 lb / 3-150 kg (11 lb / 5 kg to turn on scale) |
| Batteries | 3*1.5V AAA |
| Division | 0.2 lb / 0.1 kg |
| Unit | lb or kg |

4

4. QUESTION AND ANSWER

Why doesn't the scale work? Why doesn't the display show data? Why does the data on the screen disappear in a flash?

- Please check whether the insulation piece is taken out, batteries are installed or batteries are in low and need replacing.

Display indicates "Err" when measuring body fat.

- Your feet or scale are too wet.
- Overload the maximum capacity(330 lb / 150 kg).
- Body parameter settings should be with Height (3'3"-7'3" / 100-220 cm), Aged (10-99 years old).

The scale can not connect with App

- The app is required to be opened when weighing.
- Make sure to turn on the Bluetooth function on your device.
- Please re-turn on the Bluetooth and restart VYVO Smart App on your device.

No body fat data measured when weighing.

- Body fat scale requires weighing with dry and bare feet.
- Ensure the Bluetooth function is open.
- Please input your personal information with Height (3'3 "-7'3" / 100-220 cm), Aged (10-99 years old).


5. IMPORTANT SAFETY INFORMATION

NOTE

1. Not recommended for any person using a pacemaker or other internal device.
2. Not recommended for women who are pregnant, the measurement will be incorrect.
3. Make sure your feet are bare and dry when you measure body fat and other data.
4. Do not take measurements immediately after intense exercise.
5. Take measurements at the same or similar time of the day.
6. You can change the unit by "Setting > My information > Unit".
7. There is a maximum of 10 users per app.

CAUTION

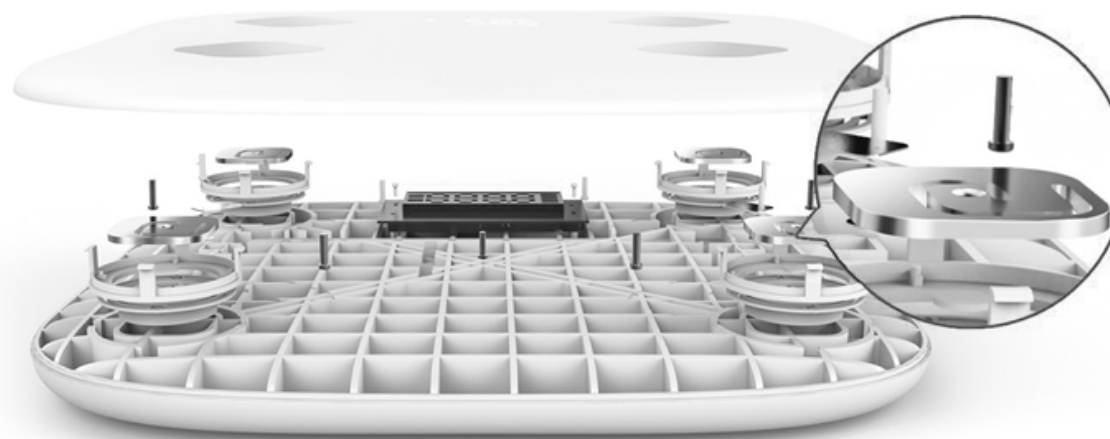1. This scale should not be used to diagnose or treat any medical condition. You should consult a your primary care doctor or physician before undertaking your diet or exercise program.
2. DO NOT stand on the edge of the scale or jump on it.
3. DO NOT overload the scale (Maximum 330 lb / 150 kg).
4. DO NOT immerse the scale in water.
5. Please use a damp cloth or a glass cleaner when you clear the surface of scale.
6. Always place scale on a hard, dry and flat surface before measurement.
7. Connect the batteries in accordance with the correct polarities indicated.
8. If you are not using the scale, it is recommended that you remove the batteries.

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